取样阀系列 Sampling valves





ENLARGE ENGINEERING CO,.LTD 134/65 Moo 2, Ban Chang Subdistrict, Ban Chang District, Rayong Province ,21130 Thailand 国家高新技术企业 浙江省著名商标 www.cndonjoy.com







目录技术资料仅供设计参考。 设计更新和技术升级恕不另行通知!

不断革新的技术、日新月异的市场发展和客户需求 是当今我们面临的最大挑战。面对这些挑战我们信 心十足。因为东正高度专业化的产品和技术创新。 是我们未来将以更高、更优的产品来满足不断增长 的市场需求。

Technological innovation, rapid development of the market and customer needs are the greatest challenges facing us today. Facing with these challenges we are confident, because of highly professional products and technological innovation, we will meet the growing demand of market and customer with our better and better products.

法律声明:

目录样册版本属东正所有。如有抄袭东正将追究 其法律责任和经济赔偿。

东正众多技术设计,都具有专利保护、加有炒菜。 东正必将追克法律责任和经济预衡。

领先源于创新 创新源于专业 东正的专业与品质的邂逅源于创新的精神



现代化工厂已经投产,期待您的光临 温州空港经济园区生产基地

Wenzhou airport economic zone production ba

占地面积29000m°,总投资1.2亿元、新引进ABB工业机器人、日本马扎克5轴加工中心、ABB机器人热光设备。

Covering an area of 29000 m2, with total investment after construction is expected to 120 million yuan, introduces new ABB industrial robot, Japan mazar-e-sharif grams of five axis machining center, ABB machine Polishing equipment.

年,注册资金5508万,是一家专注于研发高 程控制器、阀门定位器、流量控制器。阀门控制单 洁净过程控制阀、智能自控阀,高洁净泵、高洁净管 例和曾道配件。公司高度专业化的产品和技术专注于生物制 食品饮料、酿酒、精细化工等工业领域。东正众多的创 F解决方案,对满足全球,尤其是新兴市场不断增长生物制 药和食品饮料行业的需求起到了积极作用。

东正制造的产品母材料、加工技术、焊接技术、表面处理等 技术指标均已达到国内外先进水平。且符合ASME BPE、 EHEDG、FDA、3A等一系列国际标准和规范要求,先后通过 了由德国TÜV机构实施的欧盟承压设备指令(PED-97/23/EC)和(MD-06/42/EC)认证、美国 3A机构实施的3-A卫生标准认证、中国安全阀特种制造设备(TS)许可证和中国 质量认证中心CQC ISO9001:2008 质量管理体系认证等多项 国际权威认证。

东正凭借先进的技术、优异的品质和完善的售后服务体系, 先后被评定为国家高新技术企业, 省著名商标、省科技型企 业、创新型企业、专利示范企业、研发中心企业, 守合同重 信用AAA企业等众多荣誉。



关于东正 About DONJOY

Donjoy was built up in 1993, with registered capital RMB 55,08 millions. We focus on developing high quality PID process controller, the valve positioner, the flow controller, the valve control unit, high purity process control valve, intelligent selfcontrol valve, high purity pump, high purity tubing fitting and pipeline accessories. Donjoy's products and technologies are widely used in biopharmaceutical, food, drink, wine fine chemicals and other industrial fields. Donjoy's innovation solutions meet the growing demand of the industry of biopharmaceutical, food and beverage from the whole world, especially the emerging market.

Raw material, processing technology, welding technology, surface treatment from Donjoy have achieved advanced levels of domestic and foreign marekt. Donjoy productions according to ASME BPE, EHEDG, FDA and 3A standard, get certificate of PED-97/23/EC and MD-06/42/EC issued by TUV, 3A issued by 3A organization, TS licence and ISO9001:2008.

Donjoy with advanced technology, excellent quality and perfect after-sales service system, has been assessed as the national high-tech enterprises, the Provincial Famous Trademark, the Trovincial Science and Ttechnology Enterprises, Innovative Enterprises, Patent Model Enterprises, R & D Center, AAA Credi Enterprises.



创立于2005年现代化工 温州滨海经济开发区生产 velopment zone production base

DONJOY high-putity valve, pump, control system, pipe and accessories, provide the most excellent solution for various business scopes.



制造与检验 Manufacturing and testing



祝器人生产线 Robot production line

DONJOY以科学技术是第一生产力, 也是企业核心竞争力,只有不断投入核心技术研发和创新, 不断追求卓越,精益求精、以先进生产设备, 严格的质量管理体系,精密的生产检测流程, 保证每一款产品优良的性能。

全自动弯管设备 Automatic pipe bending equipment

Science and technology constitute a primary productive force, also the core competitive strength of the enterprise. Continuously intensifying the research, development and innovation of core technology, pursuing prominence and keeping improving, ensuring the excellent quality and performance of every single product by advanced production equipment, strict quality control and precise production and inspection procedures.



荣誉与资质 Honor and qualification



3A

东正是美国3A认证组织的金牌赞助商,3A是一个自愿组织,致力于维护公共卫生和如制品产品质量通过独特的行业监管程序的设备,用于处理乳制品的卫生标准3A认证。 东正通过认证产品,阀门系列55-02-1783、64-00-1784、84-02-1786、85-03-1787 53-06-1782、54-02-1580、68-00-1785,管件(63-03)证书号(1581) 等多种规格型号。



PED/97/23/EC & MD/06/42EC

东正通过了TUV 欧共体压力设备97/23/EC质量 保证模式H&材料(管件)制造认证。 阀门证书 NO.CE 0036 泵类证书 NO.705201402401-00



FDA177.2600

东正生产的阀门和泵以及流体设备等,复物制接触密封圈全部通过 或符合美国食品和药品管理局和联邦法规第24条例,CRF 21 FDA177.2600和 FDA177.1550认证,其涵盖: EPDM VITON(FPM) SILICON ANBRETFE 橡胶塑料材料等。



ISO9001/2008

东正通过中国质量认证中心CQC,ISO9001-2008/GB/T19001-2008企业 质量管理体系认证。其证书编号: 00108Q24215R1M/3302



卫生许可证

东正通过了省卫生厅卫生阀门、卫生管件等检验, 且获得了其颁发的证书,其编号: So083。 其许可证为国家出台的生活饮用水卫生监督管理方法证书。



高新技术企业

东正通过多年发展和技术创新已获得了多项国家发明专利,100多项实用 新型专利,现已成为国家高新技术企业。



浙江省著名商标

东正09年获得了温州市著名商标,2012年获得了浙江省著名商标。

DONJOY不但拥有以上认证和荣誉,且被评定为省科技型、创新型、专利示范、研发中心、守合同守信用企业等众多的荣誉证书。



工作原理 Operating principles

- ◆该无菌取样阀为每次取样前和取样后进行杀菌。
- ◆开启阀门:取样-通过旋转手柄(可调整阀门)或通过拉动操作杆(开/关阀门)来开启手动操作的阀门。然后柄和隔膜缩回, 允许液体通过。
- ◆关闭阀门: 消毒-**通过相反的**步骤,内置弹簧关闭阀门,保持软管间的通道开启,用于消毒。如果使用蒸汽,建议在出口装一个小型的压力安全阀。
- ◆标准设计:该无菌取样阀由三个部分组成、阀体、手柄和隔膜阀。橡胶膜置于阀杆上,作为拉伸塞子。
- ◆连接形式: 螺数 罐焊接、钢管焊接、卡箍。
- ◆手柄形式。(1)旋转手柄(2)可移动钥匙。
- ◆应用领域: 该阀门用于啤酒、牛奶、饮料、果汁、制药等行业。
- The aseptic sampling valve is designed for sterilization before and after each sampling.
- Open: Sampling-The manual valve is opened by rotating the handle(adjustable valve) or by activating the lever(on/off valve). The stem and the membrane are then retracted, allowing liquid to pass.
- ◆Closed: Sterilizing-The built-in spring closes the valve and keeps the channel open between the hoses though reverse process in order to sterilize. If using by steam, a small pressure relief valve(optional) on the outlet is recommended.
- Standard Design: the aseptic sampling valve consists of three parts, a valve body, a handle and a membrane. The rubber membrane is placed on the stem of the valve and works as a stretchable plug.
- ◆Connection: Thread, pipe welding, tank welding, Clamp.
- ◆Handle: (1)Turning handle (2)removable key.
- ◆Application: The aseptic sampling valve is applied to beer, milk, beverage, fruit juice, pharmacy industries tec.

材质 Material

- ◆ 阀体: AISI 316L
- ◆Valve body: AISI 316L
- ◆手柄: AISI 304
- ♦ Handle: AISI 304
- ◆膜片: EPDM PTFE
- ◆Membrane: EPDM PTFE

技术规范


◆最大压力: 600kpa(6bar)

安装尺寸 Installation dimension

气动取样阀 Pneumatic sampling valve




特点: 取样后自动关闭阀门

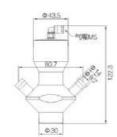

气动无菌取样阀 Pneumatic aseptic sampling valv NO.VQY100Q1

Features:

Sampling valve automati after sampling

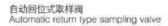



特点: 取样后自动复位



手动手轮无菌取样阀 Manual aseptic sampling valve NC.VQY100S1


Features:

Automatic reset after sampling

式取样阀 ype sampling valve



板手式手柄每次取样后可将 其取走防止意外发生

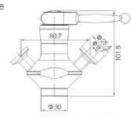

手动手柄无菌取样阀 Manual aseptic sampling valve NO.VQY100S2

This type can be removed after each sampling to avoid accident open.

温度 Temperature

- ◆最高温度: 130°C(EPDM) ◆Max. Temperature: 130°C(EPDM)
- ◆最高温度: 150°C(PTFE) ◆Max. Temperature: 150°C


Technical specification

Max. Working pressure: 600kpa(6bar)

工作原理 Operating principles



开启阀门,取样 Open the valve, sampling

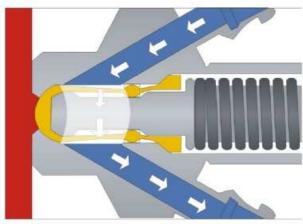

关闭阀门。消毒 Close the valve,sterilizing



无菌取样阀 Aseptic sampling valve



取样阀+管路连接器 Sample valve+pipeline connector



取样阀+罐体连接器 Sample valve+tank connector



进出口+防护套 Sample valve+protective sleeve



气动无菌取样阀 Pneumatic aseptic sampling valve



钥匙手柄无菌取样网 Lock handle sample valve



自动复位无菌取样阀 Auto reset sample valve



取样进出口选型

Interface selection:

- ◆软管Hoses
- ◆焊接Welding
- ◆卡箍Clamps
- ◆螺纹Thread
- ◆凸台口Boss Ends

罐体连接(焊接) Welded connection 罐体连接(卡箍) Clamp connection

各种配置阀体 Valve body with various configurations



密封圈 Gasket ◆ EPDM ◆ PTEE



流量微调阀头 Flow regulating valve Pneumatic valve head Auto reset valve head head



气动阀头



自动复位阀头



钥匙手柄阀头 Lock lever valve head

取样阀 Sample valve

应用 Applications

- ◆主要用于罐设备或管道里面的物料样品采集,也可安装在罐设备或管道的最低点用来排泄残留液。
- its mainly applied to extract sample from tank or pipeline. It also can be assembled at the bottom of the tank or pipeline to drain away the remains.

操作原理 Operating principle

取样阀 Sample Valve

- ◆手轮顺时针旋转阀杆边同密封圈抬起,此时阀门被开启取样, 手轮逆时针旋转阀杆连同密封圈下降此时阀门被关闭。
- •when turn the hand wheel clockwise, the valve shaft and seal rise up. At this time, valve is open for sampling. When turn the hand wheel anti-clockwise, the valve shaft and seal go down, at that time, the valve is closed.




取样阀(椭圆式手轮) Sampling valve with elliptic type handle

技术参数 Technical specification

- ◆接口方式:螺纹、卡箍、焊接
- ◆最大压力: 10bar 145psi
- ◆最小压力:真空
- ◆适用温度: -10°C/14° F 至120° F/248° F
- ◆ 阀体材质: AISI304/316L
- ◆密封材质: PTFE
- ◆Connections:thread, clamp, w
- ◆Max. Pressure:10bar 145ps
- ◆Min. Pressure: vacuum
- ◆Temperature:-10°C/14° F to 120° F/248° F
- ◆Body material: AISI304/816L
- ◆ Seal material: P7FE
- ◆O型图: NBR
- ◆ 資杆: AISI304/316L
- ◆卡簧: AISI304
- O ring: NBR
- ◆Valve shaft: AISI 304/316L
- ◆Spring: AISI304


各种连接取样阀



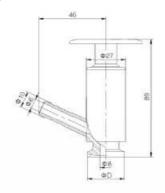

卡箍取样阀 Clamp sampling valve NO VQY210S1

SIZE	D
1"	25.4
2	50.5

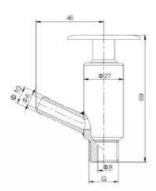

外螺取样阀 Male sampling valve NO.VQY230S1


SIZE	
1/4*	1/4°G-19
3/8"	3/8°G-19
1/2"	1/2°G-14
3/4"	3/4°G-14
1*	1'G-11



取样阀 Sample valve

直通式取样阀 Straight sampling valve



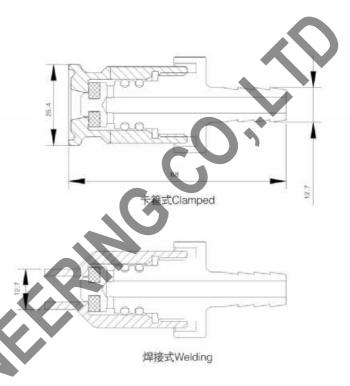



酸奶取样阀 Yogurt sampling valve

- ◆ 阀座开到罐内部或其他的设备内部,由此避免了因为管道过压 造成阀门意外开启或泄漏等情况。 ◆该阀门为常闭型(通过弹簧关闭),用
- 手板阀门开启, 松 开手柄门将自行关闭。
- ◆完全排水(液体)设计,防止了物料滞留阀门内部。
- side of tank or other equipment that it acident due to high pressure in The valve seat opens inside of can avoid opening by ac
- pipeline.

 The valve is normally closed(by spring). The valve is open as soon as pushing the handle and closed when loosing up.

 The design of completely drainage water is to avoid remains
- inside.

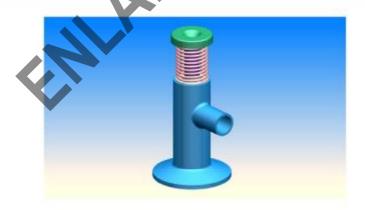

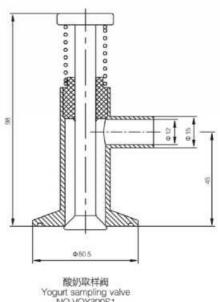

Yogurt sampling valve NO.VQY300S1




应用 Applications


- ◆该液位计主要应用在罐设备上,两个不同或相同的液位计安装 在罐设备上,可观察罐内液位体情况或观察罐内容量,也可以 直接取样。
- ◆规格: 1/4", 3/8", 1/2", 3/4"
- ◆连接:螺纹、卡箍、焊接
- ◆认证: PED 97/23/EC, FDA 177.2600
- ◆ The level gauge is mainly applied to the tank equipments. The two same gauges or two different gauges are assembled on the tank so that the liquid condition or capacity can be odserved and the sample can be extracted.
- ◆Size:1/4", 3/8", 1/2", 3/4"
- ◆Connections: thread, clamp, weld
- ◆Certificate: PED 97/23/EC,FDA 177.2600

技术参数 Technical specification


- ◆最大压力: 10bar 145psi
- ◆最小压力:真空
- ◆适用温度: -10° F/14° F至120° F/248°
- ◆ 阀体材质: AISI 304/316L
- ◆密封材质: PTFE
- ◆玻璃筒规格: Φ20
- ◆Max.pressure: 10bar 1450
- ◆Min. Pressure: vacuum
- ◆Temperature: -10° F/14° F to120° F/248° F
- ◆body material AISiS04/316L
- ◆Seal material PTF

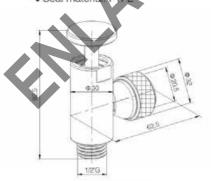

闷头液位计取样阀 Level Gauge(Superior) NO.VQY430S1-01

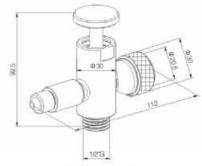

可视液位计取样阀 Level Gauge With Sample Valve NO.VQY430S1

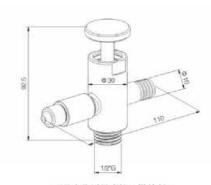

可视液位计取样阀(带软管) Level Gauge With Sample Valve(With hose adaptor) NO.VQY430S1-02















CHINA DONJOY TECHNOLOGY CO.,LTD.

DONJOY looking forward to your visit! Your telecommunications!





地址 温州龙湾区空港新区港富路25号 总机: 0577-86990888

传真: 0577-86998000

Add: No.25 Gang fu Road, kong gang new District, Longwan Wenzhou, China

Tel: +86-577-86990888 Fax: +86-577-86998000

E-mail: donjoy@cndonjoy.com

www.cndonjoy.com